we° TEMP BAND

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PWB-600T

Non contact infrared thermometer

Before using your thermometer, you should read and understand all instructions and follow all warnings.

partron

Copyright © 2020 All rights reserved.

Written on : 2020 - 10

 **Bluetooth®**

There is a possibility of radio confusion during operation of the corresponding wireless system.

/ Product introduction /



purpose of use

This product is a portable measuring device that can measure body temperature and manage information by collecting infrared radiation energy generated from the human body (forehead, under the ears, etc.) through a sensor. It is used for measurement after being located

***APP**

- After installing the 'URBAN FIT' app of Partron in the Android Play Store and iOS App Store, you can check body temperature information and other exercise information by linking with the app.
- Use an approved app downloaded from the App Store or Play Store. (For iOS, do not use the device in a randomly modified environment.)

/ Product Composition /

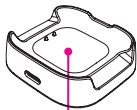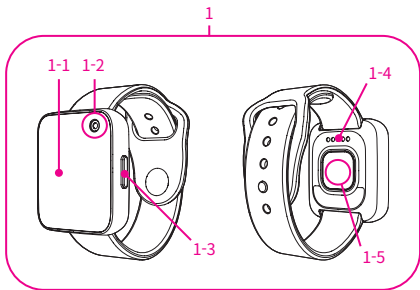2

3

4

1. Body

1-1. TFT display(1.3inch)

1-2. Temperature sensing unit

1-3. Side button

1-4. Charging terminal

1-5. Heart rate sensor

2. Charging Cradle


3. Charging cable

4. User Manual


/ How to charge the battery /

- Attach the band body to the charging cradle according to the direction. (If the direction of the charging cradle is not correct, charging will not work.)
- Appropriate KC certified power adapter must be used. (DC5V / 0.2A)


STEP 1

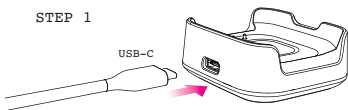STEP 2

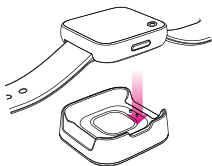STEP 3

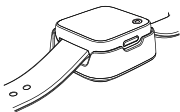/ How to wear WE-TEMP BAND /

- The product should be placed on the elbow side, avoiding the wrist bone, and fitted snugly on the wrist to measure the heart rate accurately.

/ Turning the power on and off /

STEP 1

1 click

STEP 2

- Click the side button on the side of the product once to turn on the power. When the power is turned on, it switches to the clock face.

(Time and date information is automatically updated only when connected to the APP.)

* If not linked, the clock and date will be reset.

- If there is no operation for 5 seconds, the power is automatically turned off. (In body temperature mode)
- Press the side button of the product for more than 3 seconds -> Display the power off pop-up -> Touch "OK"

/ How to operate the screen /

- Screen top, bottom, left, right slide (You can freely slide the screen up, down, left and right to select the desired function)

- Screen order

(Watch ▶ Body temperature ▶ Heart rate ▶ Workout ▶ Breath ▶ Music Message ▶ Record ▶ Setting)

/ How to measure body temperature / (The body temperature measurement function is a medical device function.)

STEP 1

STEP 2


- When you slide the screen to the left on the watch face, the body temperature measurement screen appears.

STEP 3

STEP 4

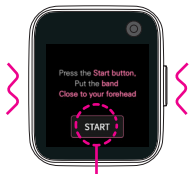Touch the start button

- After touching the body temperature measurement screen, touch the 'Start Measurement' button to start body temperature measurement.
- If the vibration sounds after touching the 'Start Measurement' button, the thermometer must be placed in the center of the forehead above the eyebrows at a distance of 3cm so that it does not touch the skin for accurate measurement.

STEP 5

STEP 6

STEP 7

- If you press the 'Measurement' button, Ready, 3, 2, and 1 are displayed in sequence, and the measurement starts.
- When the measurement is completed, a vibration sounds and the measured body temperature is displayed on the screen.

*If the measured temperature is higher than 42°C, it is displayed as HIGH, and if it is below 32°C, it is displayed as LOW.


42°C ↑

32°C ↓

- Press the 'OK' button to return to the body temperature measurement screen.

*Check the temperature measurement location and start the measurement.

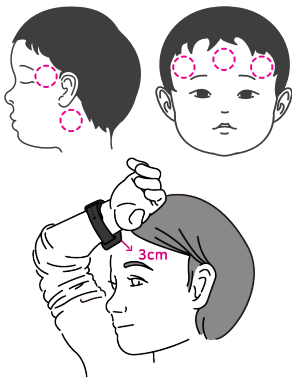- The thermometer measures the center of the forehead or the temple side with a distance of less than 3cm so that it does not touch the skin.
- If the temple or the top of the eyebrows is covered with hair, or if it is unavoidable to measure it, measure the back of the ear.
- The temperature of the left and right temples may be different. It is recommended to always measure with the higher side according to the measurement method in the user manual.

/ Changing the temperature unit /

Celsius

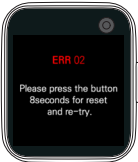
Fahrenheit

- If you press the button in the body temperature unit change part, it changes from Celsius to Fahrenheit.

/ What to do when an abnormality occurs /

	ERR 01	ERR 02
		
condition	This message is displayed when the room temperature is outside the operating temperature range of 15°C to 40°C.	This message is displayed when any mechanical error occurs.
Resolution	<ul style="list-style-type: none">- Wait for at least 30 minutes before measuring the thermometer in an available ambient temperature environment.- If the same indication appears afterwards, ask for A/S.	<ul style="list-style-type: none">- Press the side button for 8 seconds to reset the device and then try to measure again.-If the same indication appears afterwards, request A/S.

Before use,
check the following before use.


Excluding the body temperature measurement part of the WE-TEMP BAND (PWB-600T) product, the other smart band functions (heart rate measurement, exercise check, etc.) are not medical device functions.

Broadcasting Communication Equipment User's Guide

Since there is a possibility of radio wave interference during operation of the radio equipment, services related to human safety are not available.

/ Watch display /

- The watch face of the band can be changed in 4 styles, and the watch On the screen, you can check the battery / heart rate / step count / weather / BT connection status. (Different screen composition by clock type)

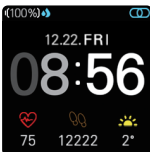(100%)

battery

Water lock

Bluetooth connection

Bluetooth disconnection

/ Snap function setting /

- Sliding down the top of the watch screen brings up the snap screen where you can set the function of the band.
- You can turn on/off notification reception / power saving mode / stretching notification / water lock on the snap screen.

/ Heart rate measurement /

- Heart rate measurement starts when you touch the screen on the heart rate measurement screen. Heart rate measurement is measured for about 30 to 40 seconds, and when the measurement is completed, the heart rate is displayed along with vibration.
- If the band is loosely worn or there is movement, heart rate measurement is not possible and a re-measurement pop-up appears. (Different screen composition by clock type)


/ Workout /

Workout

info

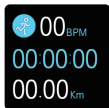- You can select an exercise mode on the exercise screen.
(Running/Walking/Hiking/Jumping Rope/Squat/Lunge/Plank/Elliptical)
- When you select the exercise mode, the countdown starts and the exercise information screen appears.
- If you press the end exercise button at the bottom of the exercise information screen, the exercise ends.

/ Auto workout mode /

- In the case of walking and running, it automatically enters the exercise mode by detecting movement for a certain period of time without any additional settings. (If the number of steps detected is 100 steps for 1 minute for women (5 minutes maintenance), 95 steps for 1 minute for men (5 minutes maintenance))
- When exercising by selecting the exercise mode while running the automatic exercise mode the automatic exercise mode ends.

Snap screen


Show icon on launch

- * When the automatic exercise mode is detected and measured, the accumulated data on the main screen of the band can be checked in the APP (check in the past exercise history view after the exercise ends)

/ Breath /

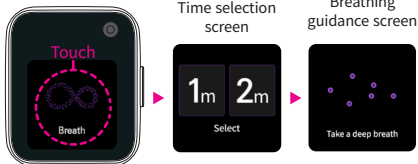- You can select the breathing time by touching the screen. It vibrates every 5 seconds to guide your breathing, and you can take a deep breath by following the guide.

/ Message /

- On the message screen, you can check the received message by touching the screen.
- The messages are displayed in the order of receiving SNS/ SMS / APP PUSH notifications.
- Only up to 5 received messages are displayed.
(The old messages are automatically deleted, and the latest messages are saved.)

/ Music /

- When you touch the screen on the music screen, you can play music from the terminal connected to the band.
- Music playback can also be used in exercise mode and cannot be used if there is no terminal connected to the band.
- After connecting the band to the APP, the music player must be activated when the music is played for the first time in order to play the music normally. (Run APP → Click Music → Activate Music Player → select the desired media app)

/ Play music /


- Music playback is also available in exercise mode, and you can play music during exercise by tapping Start Music at the bottom of the exercise mode.
- All music playback functions cannot be used if there is no terminal connected to the band.

/ Settings /

- If you touch the screen in the setting screen, the detailed setting screen appears.
- On the detailed setting screen. You can set (Change Watch Face / Adjust Screen Brightness / Raise Wrist / Reset Band).

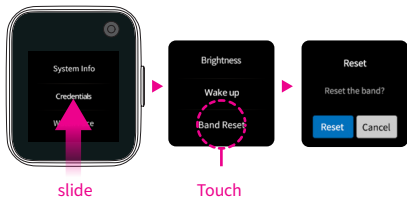- If you slide the screen from the bottom up on the detailed setting screen, the add detailed setting screen appears.
- You can reset the band by pressing the band reset button on the detailed setting screen.

/ APP connection /

(To connect with a smartphone, you must use a dedicated APP.)

* When using it for the first time, please use it after connecting it to the Urban Fit APP.


1. Search and install Urban Fit in the Android Play store or iOS APP store.
2. After running the APP, agree to the terms of use and set up a profile.
3. When you have finished setting up your profile, tap Find Band in APP.
4. When the band is found, the PIN number will appear on the Urban HR band's screen.
5. Enter the PIN number of the band in the APP and press OK to connect to the APP.

/ APP function /

(You must use a dedicated APP for connection with a smartphone.)

1. Record check



When running the APP, it is connected to the record check mode.

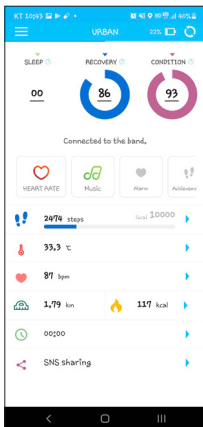After the band and the APP are connected, the recording of heart rate, exercise, body temperature, etc. can be checked in the APP. If the band and APP are not connected, the APP cannot check the record.

2. Workout mode

Scrolling left in record check mode moves to exercise mode.

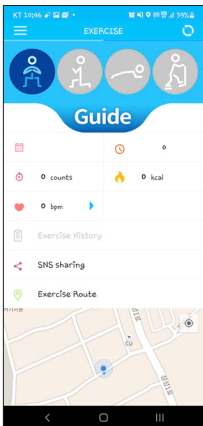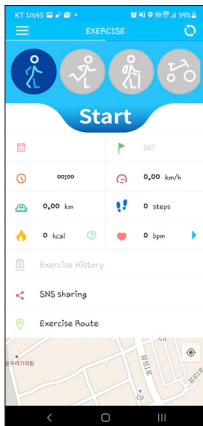
With the band and the app connected, you can do walking, running, hiking, and cycling workouts, and you'll get a workout posture guide for jump rope, squat, lunge, plank, elliptical workouts.

* In the case of cycling, you can check your running and exercise records only in the app.

3. Settings

When selecting the upper left (☰) menu, it moves to the main menu screen.

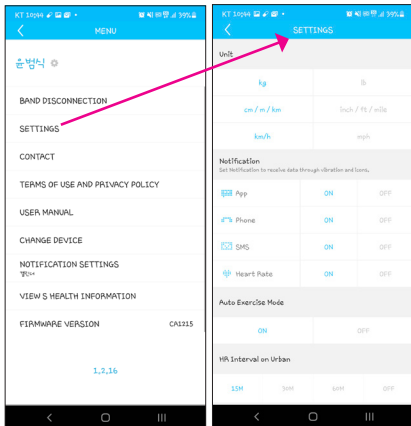- ① Unit setting: It is possible to set the unit of exercise record.
- ② Notification settings: You can set the notifications you want to receive phone calls, text messages, social media, schedule notifications, and heart rate. "For iOS, App notification -> SNS (KakaoTalk) notification provided"
- ③ Automatic exercise mode setting: You can set the automatic exercise mode.
- ④ Cycle setting: You can set the normal heart rate measurement and the heart rate measurement period during exercise.

Precautions for safety

1. Do not disassemble or modify the main body. It may cause fire or breakdown.
2. Do not use or store in humid places, places with volatile substances, or near fire.
3. If you use a car charger, it can be used up to a general vehicle.
4. If multiple apps are used at the same time, the app may be automatically closed by the mobile operating system, and in this case, the Bluetooth connection may be disconnected.
5. 'Urban fit' application gets relevant information online to measure distance accurately. Data charges may apply when connecting via 3G or 4G during this process.
6. The heart rate measurement function of this product is used only for physical fitness and health, and cannot be used for diagnosis of diseases or other conditions, treatment or prevention of diseases.
7. This product supports a certain level of water resistance, but does not guarantee complete water resistance under special conditions (swimming, diving, diving). Also, manufacturers and distributors cannot be held liable for errors arising in these circumstances.
8. Be careful not to let the terminals of the charging cradle connected to the charging cable come into contact with metal objects. It may cause fire or malfunction.

/ Product management method /

1. Do not use and store at high or low temperatures.
2. Avoid fire and direct sunlight, and avoid humid places or pollutants. Please keep it out of place.
3. For colored cosmetics (ex, female color cosmetics such as sunscreen) If it is contaminated, you should wipe it off immediately with wet wipes. If the contamination is left for a long time (30 minutes) or more, the stain may not be erased, so use with caution.
4. Water resistance and performance may be deteriorated if the product is dropped or impacted.
5. Use a heating device such as a hair dryer to clean the product. Do not dry.
6. Since the temperature sensor part is a very sensitive part, accurate measurement Please keep it clean and undamaged at all times.
7. When the sensor and the body are contaminated, a 70% diluted solution of ethyl alcohol is applied. Carefully wipe with a cotton swab or alcohol swab. (Recommended to use after about 10 minutes after cleaning the sensor)

/ Maintenance and management /

1. When using multiple apps or functions at the same time, the battery may drain quickly, and the power may be turned off during data transmission, so always fully charge the battery before using.
2. If the power does not turn on or does not work even after charging the band, contact the customer service center (1644-4068) or ask for A/S.

/ Quality assurance standards /

1. The warranty period of this product is one year from the date of purchase.
2. Even within the free warranty period, if the service is eligible for a paid service, a paid fee may be charged.
3. Only the purchaser is guaranteed the product, and you can request proof of it.
4. If the product label is removed, physical changes to the product, or repaired by an unauthorized third party, the warranty will be void.
5. Product specifications and appearance are subject to change without prior notice for quality improvement.
6. The cleaning of the charging terminal must be managed periodically, and if corrosion or perforation occurs due to negligence, exchange or return is not possible.
7. The smart band supports a certain level of water resistance (IPX8), but in the case of damage to the product due to water immersion, a fee may be incurred.

Waterproof level guide and precautions for use

This product supports water resistance in life, but pay attention to the following for normal waterproof operation.

- Be careful not to let the product come into contact with high water pressure.
- Do not use when diving, diving, snorkeling, flow and pressure water sports, etc.
- If the device or hands are wet, dry them before use.
- If the product is exposed to liquids other than clean water (salt water, soapy water, oil, perfume, detergent, cosmetics, etc.), wash it with clean water under weak water pressure and wipe it with a soft cloth before use.

* This product has been tested in a controlled environment and It supports waterproofing under conditions.

Type of occurrence	Provided details	
	Within the warranty period	After the warranty period
1. Within 1 month of purchase Defects in performance and functionality occurred during normal use	Product exchange or Warranty or Refund (purchase site)	Not applicable
2. Partron Service If the same symptom occurs more than 3 times after repair by an engineer	Product exchange or Free repair	Paid repair
3. User inexperienced And failure due to carelessness (flooding, shock, damage, unreasonable operation, etc.)	Paid repair	Paid repair
4. In case of failure due to user's arbitrary modification, disassembly or repair	Paid repair	Paid repair
5. In case of failure due to the use of consumables and options not provided by Partron (protection film, strap, etc.)	Paid repair	Paid repair
6. In case of failure due to natural disasters (fire, earthquake, flood damage, tsunami, etc.)	Paid repair	Paid repair

※ We do not offer compensation for product failure or damage.

/ Product disposal /

This product is a product that does not generate any other harmful residues to the human body, and in principle, it is separated and discharged according to general electronic product disposal procedures.

/ Transport and storage conditions /

- Temperature : -25°C ~ 55°C
- Humidity : 15% RH ~ 90% RH
- atmospheric pressure : 70.0 KPa ~ 106.0 KPa

/ Operating condition /

- Temperature : 15°C ~ 40°C
- Humidity : 15% RH ~ 90% RH
- atmospheric pressure : 70.0 KPa ~ 106.0 KPa

Refer to instruction manual


Keep dry


Keep away from sunlight

Temperature limit

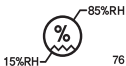Humidity limitation

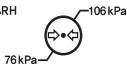Atmospheric pressure limitation

Caution

Type B Applied part

Recycling

R-C-Ptr-PWB-600T
Broadcasting communication equipment, etc.
Conformity assessment indication

Transmission data(Work out information)

/ Product specification /

Specification	WE - TEMP BAND	
Item Name	PWB-600T	
Manufacturer	Patron Co., Ltd.	
Battery	DC 3.7V 140mAh, Lithium polymer rechargeable	
Display	1.3inch TFT display	
Device Classification / Waterproof Rating	Continuous operation, IPX8	
Resolution	0.1	
Compatibility	Android 5.0 or higher / iOS 10.0 or later	
Sensor	heart rate, 6-axis acceleration	
Dimension (WxLxH)	36.23 x 43 x 12.25mm (without strap)	
Weight	36.5g (including strap)	
Power adapter rated powerw	DC5V / 0.2A	
Product life	3 years	
Voltage Rating	Type of protection against electric shock	Internal power type
	Degree of protection against rating	Type B device
Bluetooth	BLE 5.1	
Software	PWB-600T	
Manufacture Date / Lot Number	Displayed in the label on product package	

/ warranty /

Product Name	WE - TEMP
Model	PTD-600T
Manufacture Date / Lot Number	Displayed in the label on product package
Place of Purchase	
Date of Purchase	
Free warranty period	within one year from the date of purchase

This product is manufactured through our strict quality control and inspection process.

- Follow the “Consumer Damage Compensation Regulations” of the Ministry of Finance and Economy’s notice of compensation standards for product repairs and exchanges.
- If normal use fails, we will repair it free of charge during our service warranty period.
- If a problem occurs with the product during the warranty period, please provide the model name, serial number, purchase date and problem.
- In the event of product defect damage caused by misuse or neglect of product management, the manufacturer or agent shall not be held liable for any damage or damage.

partron

Manufacturer | partron Co.Ltd

Partron Address | 22, Samsung 1-ro 2-gil,
Hwaseong-si, Gyeonggi-do,
Republic of Korea (18449)